
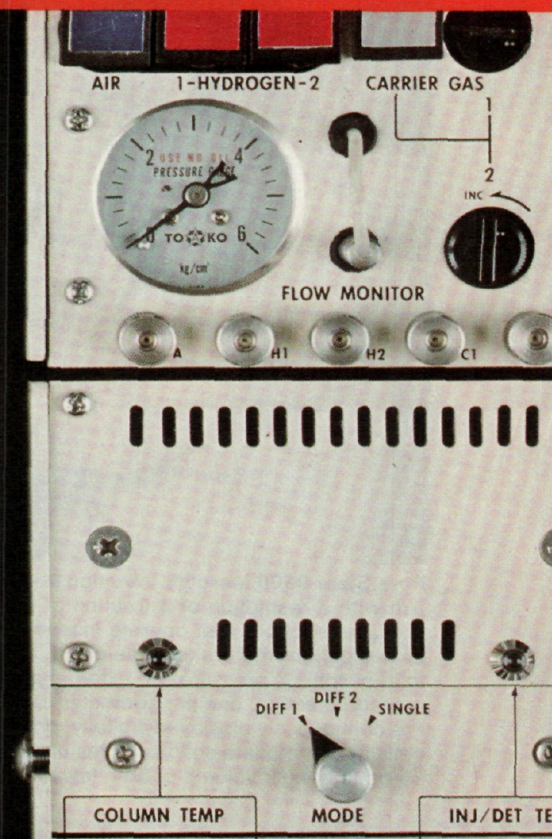


NEED THEN

**High Resolution CHROMATOGRAPHY
with capillary column?**

**Why not try the SHIMADZU GC-MINI 2
series GAS CHROMATOGRAPH?**

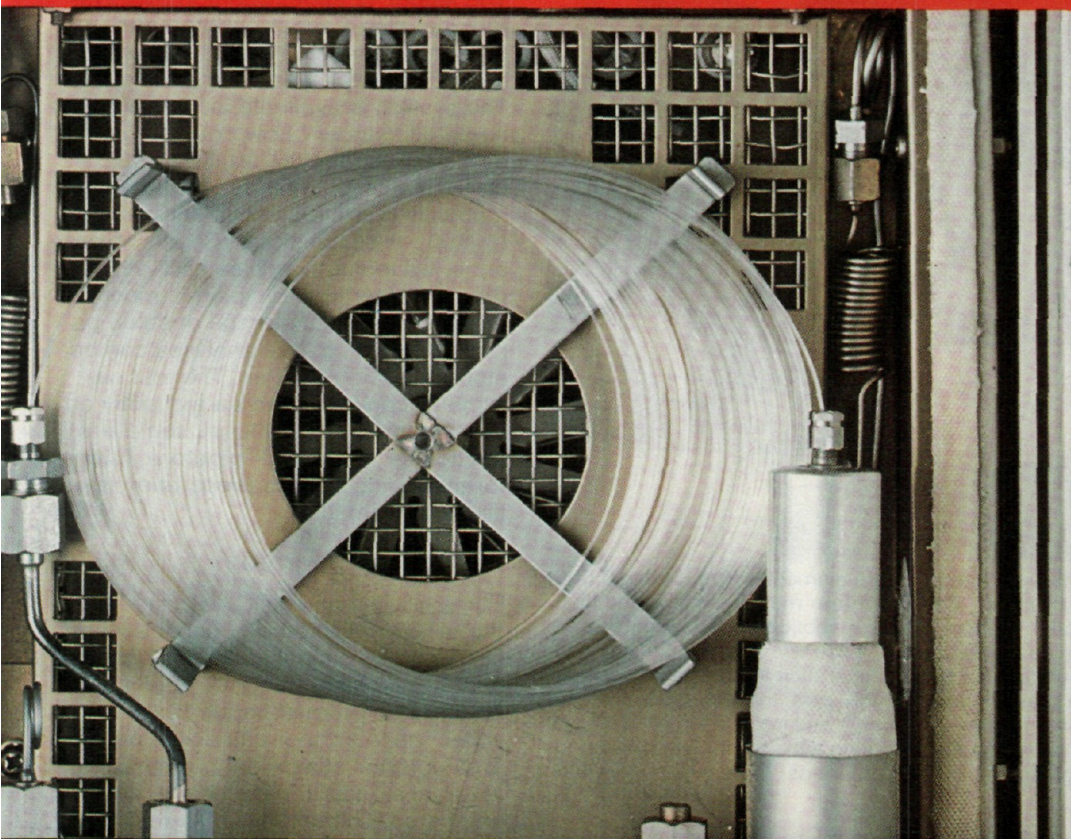

The GC-MINI 2 Series are compact, high-quality gas chromatographs specially developed for capillary column analysis.

• **Easy Column Connection**

All capillary columns, such as fused quartz and soft glass, can be installed as easily as conventional packed columns and can be tightened with common open-end wrenches. Of course, no reshaping of the column ends is required.

• **Heated Variable Splitter**

The splitter temperature is set as same as the injection port. The splitter can be maintained at the same temperature as the column oven by detaching the heat conductive mount from the splitter. The splitting ratio can be easily ad-

justed by means of a needle valve outside the oven.

• **All-Glass System**


The splitter and make-up gas system is made of glass to avoid decomposition of the sample components on the metal surface.

• **Outstanding Performance**

A closed oven system with 3 different temperature sensors ensures unbelievable accuracy and repeatability. The automatic external air intake system controls the column temperature near the ambient temperature.

A contamination-free flame ion-

ization detector and a constant current electron capture detector are available.

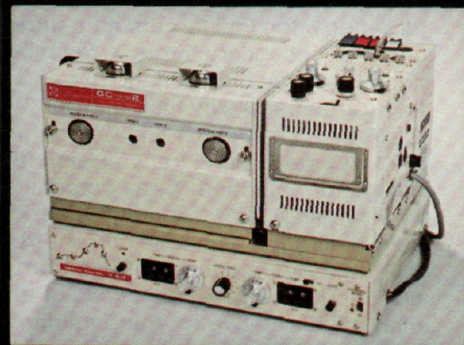

GC-MINI 2 FID model

• **Write today for more information on these and other Shimadzu instruments.**



SHIMADZU
SCIENTIFIC INSTRUMENTS, INC.

SHIMADZU SCIENTIFIC INSTRUMENTS, INC.

9147-H Red Branch Road, Columbia, Md. 21045, U.S.A. Phone: (301) 997-1227 Telex: 87959

SHIMADZU (EUROPA) GmbH

4000 Dusseldorf 1, Johannes Weyer Str. 1 F.R. Germany Phone: (0211) 314061 Telex: 08586839

SHIMADZU CORPORATION INTERNATIONAL MARKETING DIV.

Shinjuku-Mitsui Building, 1-1 Nishishinjuku 2-chome, Shinjuku-ku, Tokyo 160, Japan
Phone: Tokyo 03-346-5641 Overseas Telex: No 0232-3291 SHMDT J